

PHYS 360A/B
Experiment 20: Nuclear Magnetic Resonance

Mikhail Roslikov Kamalesh Paluru

October 17, 2023

Contents

1	Introduction	4
2	Theoretical Background	5
3	Experimental Design and Procedure	6
4	Analysis	7
4.1	Finding Resonance	7
4.1.1	Voltage at $t = 0s$	7
4.1.2	Free Induction Decay	9
4.2	The Free Induction Decay and T_2	9
4.2.1	Finding T_2	9
4.3	Measurement of T_1	11
4.3.1	Calculating T_1	12
4.4	Measurement of T_1 with Inversion-Recovery Sequence	13
5	Conclusion	14

Abstract

1 Introduction

Hello World!

2 Theoretical Background

3 Experimental Design and Procedure

4 Analysis

4.1 Finding Resonance


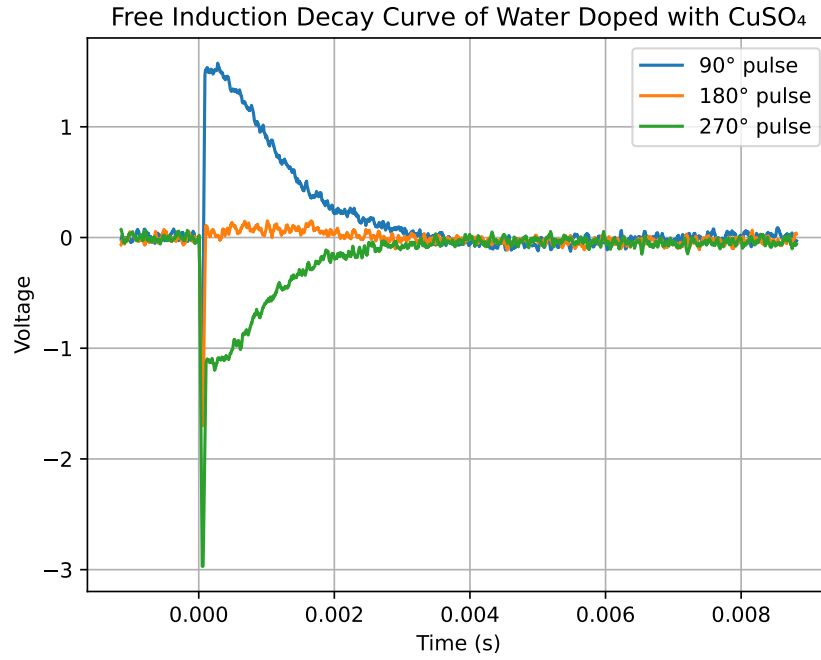


Figure 1: Free Induction Decay NMR signals for 90°, 180°, and 270° pulses

4.1.1 Voltage at $t = 0$ s

90° Pulse

- A 90° pulse is a pulse that is applied long enough to tip the magnetization vector by 90° from its initial direction (at a small angle with the positive z-axis) in the rotating frame:

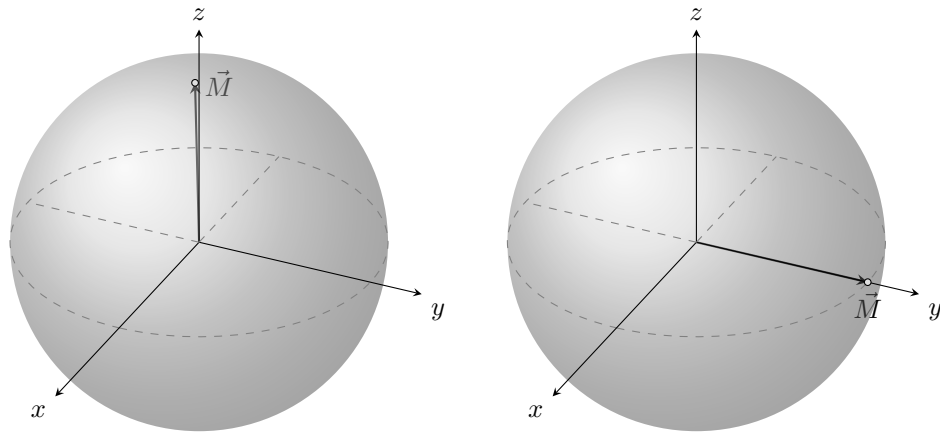


- Now, nearly half the spins are in the “up” state and the other half are in the “down” state.
- Since \vec{M} is in the x-y plane, the z component of the magnetization vector vanishes:

$$M_z = 0$$

- This is a higher energy state than the equilibrium state with the magnetization vector, \vec{M} , pointing along the positive z-axis.
- The receiver gain was adjusted so that the precessing \vec{M} induced¹ a current in the coil as shown in the 90° trace of Figure 1.

180° Pulse


- If we apply the pulse for twice as long (increase the pulse width to twice that of the 90° pulse), we rotate the magnetization vector by 180° :

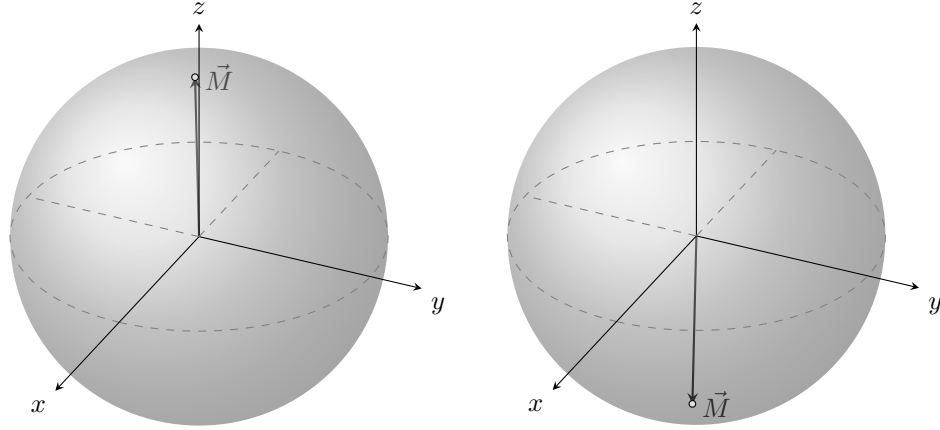

- Most of the spin are now in the “down” state.
- However, this magnetization vector doesn’t induce current in the coils since the component in the x-y plane is nearly 0.

270° Pulse

- This time, the pulse is applied long enough to rotate \vec{M} by an additional 180° from the 90° case:

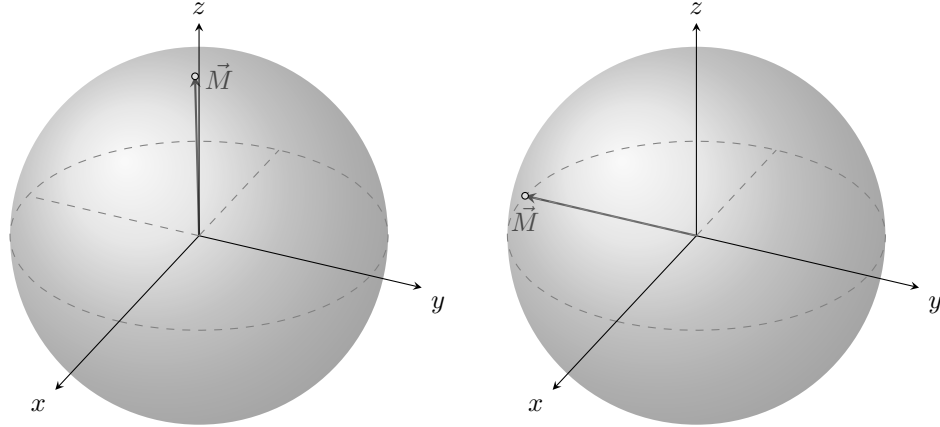

- That is, \vec{M} returns to the x-y plane but is anti-parallel to \vec{M} in the 90° pulse case:

$$\vec{M}_{270^\circ} = -\vec{M}_{90^\circ}$$

- It hence precesses in the opposite sense of rotation as the 90° case².

¹According to Faraday’s law.

²And vice versa.

- According to Faraday's law, the direction of the current \vec{M}_{180° induces in the coil is opposite to that of \vec{M}_{90° .
- This is why we see a current that is $-I_{90^\circ}$ induced in the 270° case in Figure 1.

4.1.2 Free Induction Decay

- For all three pulses, we see that the signal vanishes over time.
- Recall that for the 90° and 270° pulses, the magnetization vector is in the x-y plane.
- Because of small variations in the magnetic field that the magnetic moments, $\vec{\mu}$, for each particle experience, the magnetic moments being to randomly dephase.
- They spread out in the x-y plane causing the magnetization vector and hence the induced current to vanish as a whole.

4.2 The Free Induction Decay and T_2

- In this section, we take a closer look at the FID observed in the 90° and 270° traces of Figure 1 (for different samples).
- For the same 90° pulse, we have the following FID traces:

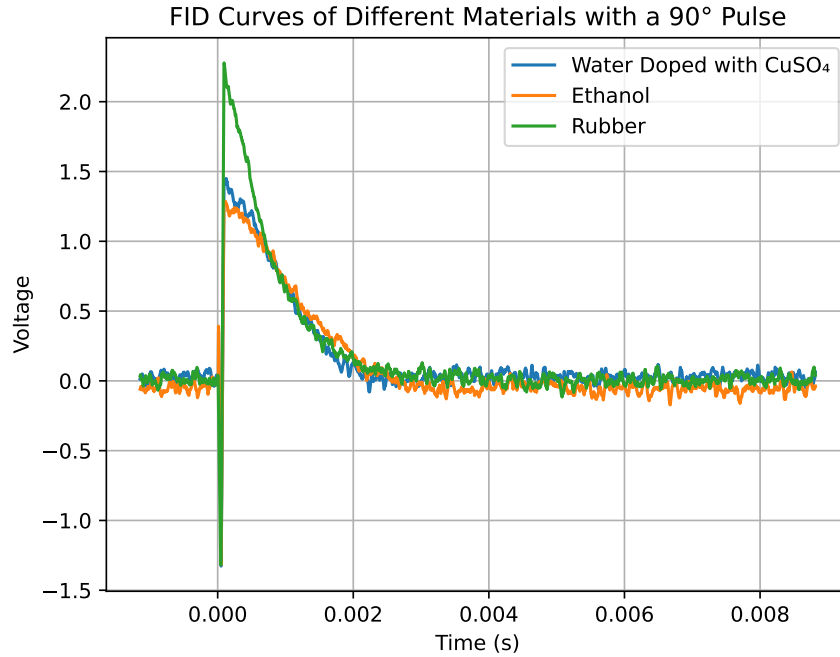


Figure 2: FID for the 90° pulse in water doped with CuSO_4 , ethanol, and rubber

4.2.1 Finding T_2

Analytically

- On plotting the traces in Figure 2 on a semi-log graph (along with the best fit lines), we get:

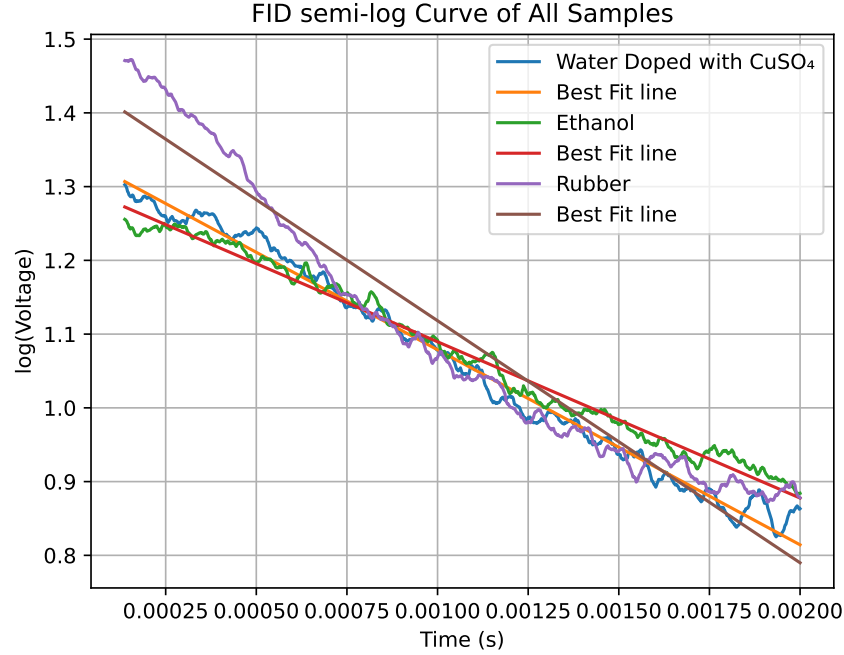


Figure 3: $\ln(M_{xy}) = -\frac{t}{T_2^*} + \ln(M_0)$ plotted for the three traces in Figure 2

- We can find T_2^* using the following equation:

$$M_{xy}(t) = M_0 e^{-\frac{t}{T_2^*}}$$

Here, $M_{xy}(t)$ is the trace seen in Figure 2,
 M_0 is the magnetization at $t = 0$,
 t is the time.

- We can take the \ln of both sides to get:

$$\ln(M_{xy}) = \ln(M_0) + \ln\left(e^{-\frac{t}{T_2^*}}\right)$$

$$\ln(M_{xy}) = -\frac{t}{T_2^*} + \ln(M_0)$$

- Comparing this to $y = mx + c$, we see that the slope, m , is given by:

$$m = -\frac{1}{T_2^*}$$

$$\Rightarrow T_2^* = -\frac{1}{m}$$

Sample Calculation

- Consider $m = -328.236\text{s}^{-1}$ for Doped Water in Table 1:

$$T_2^* = -\frac{1}{-328.236\text{s}^{-1}}$$

$$T_2^* = -\frac{1}{-328.236}$$

$$T_2^* = 0.003047\text{s}$$

	Slope	T_2 (s)
Doped Water	-328.236	0.003047
Ethanol	-264.426	0.003782
Rubber	-211.89	0.004719

Table 1: Slopes of best fit lines and calculated T_2^*

- We can hence find T_2^* from the slopes of the best fit lines in Figure 3 for the other samples:

Estimation

- We can find an estimation for T_2^* directly from the FID traces.
- For $t = T_2^*$, we have:

$$M_{xy}(T_2^*) = M_0 e^{-1}$$

$$\frac{M_{xy}(T_2^*)}{M_0} \approx 0.37$$

$$M_{xy}(T_2^*) \approx 0.37 M_0$$

- We can estimate the time at which $M_{xy}(T_2^*)$ is approximately $0.37 \cdot M_0$:

	M_0	$0.37 M_0$	T_2 (s)
Doped Water	1.46	0.540	0.00112
Ethanol	1.40	0.518	0.00122
Rubber	2.40	0.888	0.00089

Table 2: Estimate of T_2^* using the initial magnetization

4.3 Measurement of T_1

- In this section, we perform an Inversion Recovery (IR) experiment to find the spin-lattice relaxation time, T_1 .
- We use two 90° pulses $t = \tau$ apart:

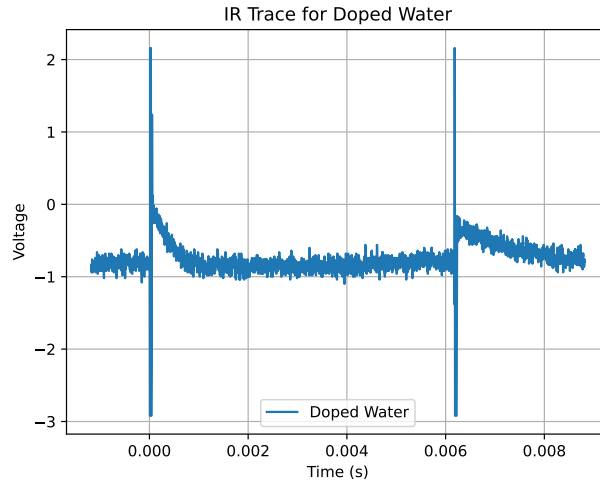


Figure 4: IR trace for water doped in CuSO_4

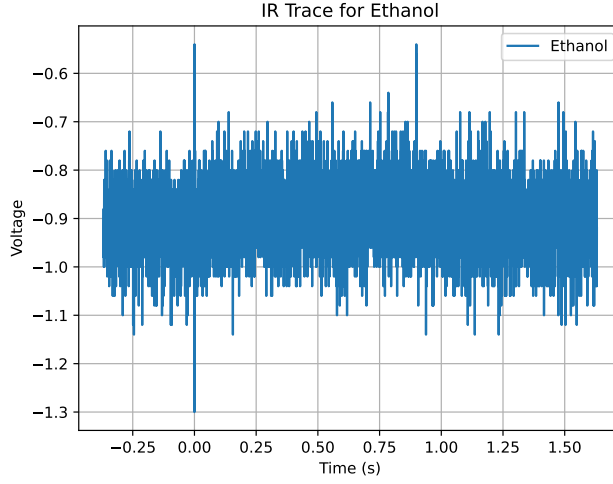


Figure 5: IR trace for Ethanol

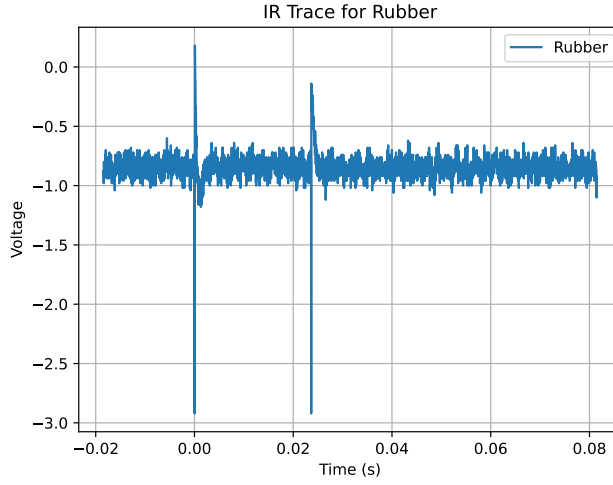


Figure 6: IR trace for Rubber

4.3.1 Calculating T_1

- To find T_1 , we use the following equation:

$$M_z(\tau) = M_0 \left(1 - e^{-\frac{\tau}{T_1}} \right)$$

- Since we already chose to record the data such that $M_z(\tau) = \frac{M_0}{2}$:

$$\frac{M_0}{2} = M_0 \left(1 - e^{-\frac{\tau}{T_1}} \right)$$

$$\frac{1}{2} = 1 - e^{-\frac{\tau}{T_1}}$$

$$\ln \left(e^{-\frac{\tau}{T_1}} \right) = \ln \left(\frac{1}{2} \right)$$

$$-\frac{\tau}{T_1} = -0.693$$

$$\frac{\tau}{T_1} = 0.693$$

$$T_1 = \frac{\tau}{0.693}$$

- Where τ is the distance between peaks³:

$$- \tau_{\text{doped H}_2\text{O}} = 0.006205\text{s} - 0.000030\text{s} = 0.006175\text{s}$$

$$- \tau_{\text{ethanol}} = 0.898870\text{s} - 0.000248\text{s} = 0.898622\text{s}$$

$$- \tau_{\text{rubber}} = 0.023670\text{s} - 0.000088\text{s} = 0.023582\text{s}$$

- And so $T_1 = \frac{\tau}{0.693}$ becomes:

	τ (s)	T_1 (s)
Doped Water	0.006175	0.008911
Ethanol	0.898622	1.296713
Rubber	0.023582	0.034029

Table 3: T_1 calculated from τ

4.4 Measurement of T_1 with Inversion-Recovery Sequence

-

³This was done by finding where the peaks happen in the data.

5 Conclusion